

OUTPERFORM - OUTSMART - OUTLAST

THE ULTIMATE HIGH SPEED TREADMILL

The ELG has everything you need for the ultimate performance lab – superior speeds, extreme elevation and endless power. Sizing at 27' wide and 96' long this running machine was built like a tank and is sure to push even the world's most elite athletes to the limit. Max out your speed at 25 mph or try some intense uphill at a max elevation of 35% or even downhill training at -5%. Push the envelope in your performance training with the diverse ELG treadmill.



ELG




- Standard Feature
- Optional Feature

		ELG
USER INTERFACE	Coded/Non-Coded Polar Heart Rate Pick-up	●
	ANT+ Heart Rate Pick-up	
	Touch Grip Heart Rate Pick-up	●
	Personal Trainer Display Board	●
	RS 232 Serial Interface	○
	CSAFE Communication Compatible	●
	Wireless Remote	○
	PT Pro Software	○
	RFID System	○
	Running Surface Dimensions	27" W x 96" L (70 x 244 cm)
SPECIFICATIONS	User Weight Capacity	Run: 400 lb. / Walk: 800 lb. (4 mph max)
	Belt Type	87 Individual Slats
	Drive System	158 Precision Ball Bearings with 16 Roller Guides (4 mm lateral tolerance)
	Running Surface	Vulcanized Rubber (38-43 shore hardness)
	Drive Motor	5 hp Continuous (15 hp peak) Brushless Servo
	Unit Weight	1,634 lb. (Shipping Weight 1850 lb.)
	Width	48" (122 cm)
	Length	102" (259 cm)
	Height	79" (201 cm)
	(-5%) - (+35%) Incline	●
	0-25 MPH (0-40 km/h)	●
	Reverse	○
	Power Supply	208/230 Vac 30 Amp Power Supply **Dedicated Circuit and NEMA L6-30R Outlet Receptacle Required
	Off-Mount Display	○
Single Handrail	○	
Gantry System	○	
Warranty	3 Year Drive, Motor, & Belt 3 Year All Components 1 Year Labor	



CONVENIENT USER CONSOLE

- Comfortable handrail control of speed, elevation, and stop
- Multiple LED readouts monitoring speed, incline, distance, time, and heart rate
- CSAFE fitness communications compatible

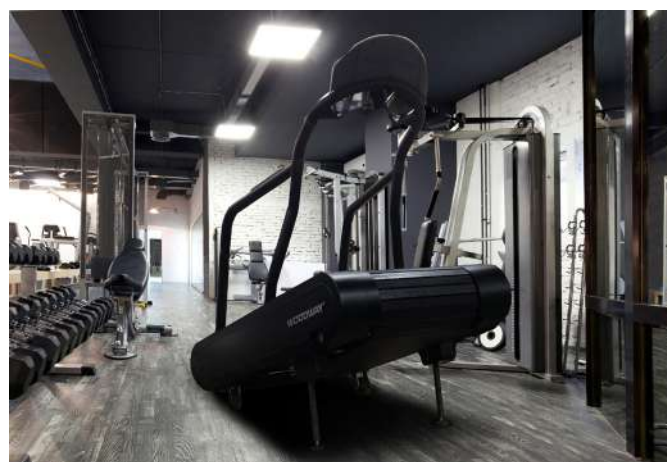

WOODWAY SLAT BELT RUNNING SURFACE

- Individual slats with vulcanized rubber surface
- Zero stretch & Zero slip toothed, wire-reinforced lateral belts
- Near frictionless precision ball bearing rail guide

SAFETY & RELIABILITY

- Emergency stop magnet and safety lanyard
- Polar monitor circuitry – integrated heart rate monitoring
- Advanced Slat Belt shock absorption
- Life expectancy of running surface 150,000 to 200,000 miles
- Extended warranty and factory renewal program

Treadmill specifications subject to change without written notice. Warranties outside the U.S. may vary - Please contact your dealer for details.



OPTIONAL SAFETY GANTRY

Reduce fall risk with an overhead gantry system that utilizes a harness to suspend individuals above the running surface in the event of a fall. The gantry includes one harness (any size) and rope suspension system.

FOR ADDITIONAL INFORMATION ABOUT THIS PRODUCT VISIT WWW.WOODWAY.COM OR CALL 800.WOODWAY TO TALK TO A REPRESENTATIVE.

WOODWAY declines any liability for injury and/or property damage which could have been avoided with use of a gantry safety system for high speed training applications.